



# One Instrument - Three Configurations

Which Analyzer is **Just Right** for Your Laboratory?



**AIA-900**  
Loader Model



**AIA-900**  
with 9 Tray Sorter



**AIA-900**  
with 19 Tray Sorter

The latest addition to Tosoh's family of AIA immunoassay analyzers is available in three configurations, allowing you to choose a system that's **just right** for your laboratory workload.

**AIA-900**  
AUTOMATED IMMUNOASSAY ANALYZER

- Instrument can be re-configured as testing volumes increase
- Throughput of 90 tests per hour
- Utilizes Tosoh's Unit Dose Test Cup reagent technology



**TOSOH BIOSCIENCE**


# A System that Grows as Your Testing Volume Increases



**AIA-900**  
Loader Model



**AIA-900**  
with 9 Tray Sorter

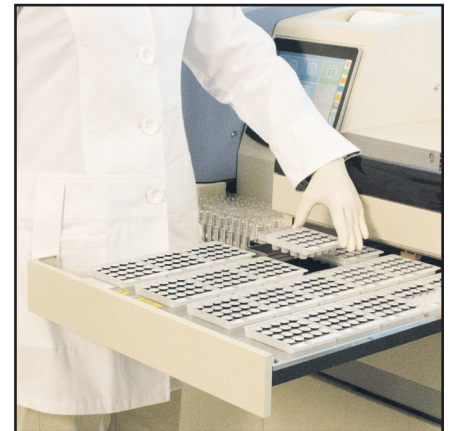


**AIA-900**  
with 19 Tray Sorter

## Unit Dose Test Cup Technology

- Single, unitized cups use a dry reagent format that requires no pre-mixing, no pre-measuring, no on-board refrigeration and no waste.
- Results are interchangeable with all Tosoh AIA Systems.
- Dry reagent format ensures calibration stability of up to 90 days for most assays.
- Reagents have a one year shelf life from date of manufacture.
- The AIA test cups and trays are labeled with code and lot number for automated scheduling and inventory management on AIA-900 sorter models.
- ST AIA-PACK reagents provide a faster turnaround time due to a 10 minute incubation protocol.
- Extensive test menu\* provides a wide range of immunoassay tests including thyroid, cardiac markers, tumor markers, reproductive hormones, anemia markers, metabolic markers, diabetes, kidney markers and more.

\*AIA-900 is pending FDA clearance for ST AIA-PACK PA

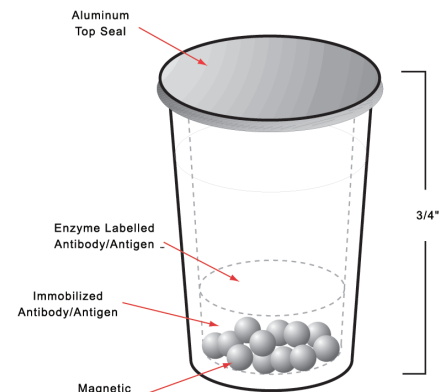

**Test Cup Diagram**



**TOSOH BIOSCIENCE**

Tosoh Bioscience, Inc.  
6000 Shoreline Court, Suite 101  
South San Francisco, CA 94080  
Tel: (800) 248-6764  
Fax: (650) 615-0415  
[www.tosohbioscience.us](http://www.tosohbioscience.us)